



Access Health

Athelgen Apollo Chiropractic Table

APOLLO CHIROPRACTIC TABLE

The Apollo 5 Advantage features Thoracic and Lumbar magnetic drops. With fewer moving parts, magnetic drops last longer and provide a smooth, crisp and responsive action. Our new drops have improved strength and have a lower minimum hold resistance.

The improved Ultra firm Australian Dunlop cushioning system provides superior patient comfort and is perfectly suited to chiropractic adjustments and manipulations. An extra powerful German motor enables fast height adjustment with less wasted time during treatments.



ATHLEGEN

No assembly
is required!

S.W.L
220 kg

KEY FEATURES

- Magnetic Thoracic and Lumbar drops
- Extendable footrest, handy for taller patients
- Height adjustable footrest for greater patient positioning flexibility
- Advantage contoured top
- A built-in sturdy adjustable head section with armrest
- Paper roll holder and cutter
- Choice of PVC, PU and genuine leather fabrics
- Choice of two Dunlop cushioning systems
- Proudly made in Ballarat, Victoria for more than 25 years

15 Year
Warranty
Structural

5 Year
Warranty
Hardware

3 Year
Warranty
Motor

Ph: 03 9878 1399

| sales@accesshealth.com.au




| www.accesshealth.com.au



Access Health

Athlegen Apollo Chiropractic Table

SPECIFICATIONS

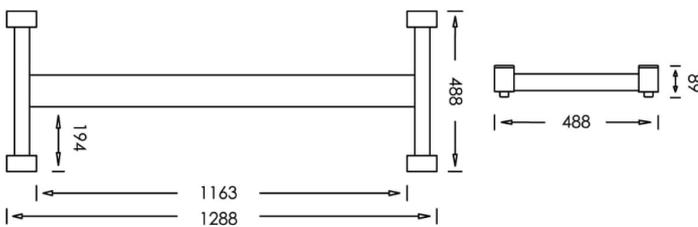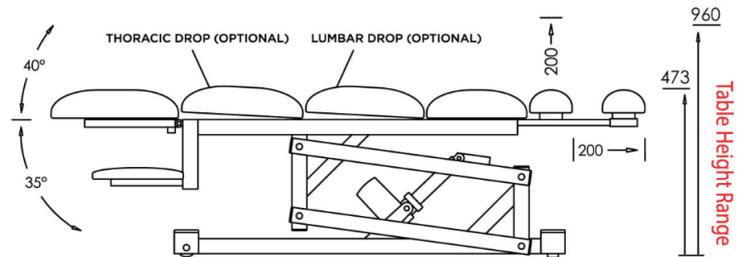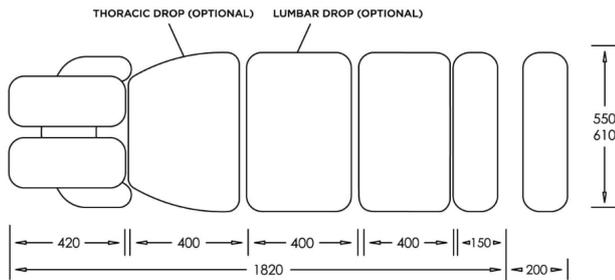

Motor	Lift force: 10,000N Speed: 37mm/sec Height adjustment (min to max): 487mm in 13 secs Speed test load: 180kg
Footswitch	Ergonomic low profile design (20mm) Stainless steel Pneumatic control
Power Requirements	240VAC @ 3.4 amps Power cord: IEC 3 pin, 10 amps (3m)
Height Range	473 to 960mm
Table Width	550mm, 610mm or Advantage Top (490-550mm)
Table Length	1820 to 2020mm

Cushioning	Dunlop® Australia cushioning Ultra-Fresh® antimicrobial treatment (CFC-free) Polyurethane foam
Upholstery Fabric	Evolve polyurethane (PU) or Encore PVC Water-based formulas (90%) Made from polycarbonate resins Low VOC No heavy metals Formaldehyde-free
Upholstery Board	PYNEboard® E0 Made from renewable plantation pine Low formaldehyde emission board
Base Style	Concealed in-line wheels
Shipping Dimensions	1930 x 730 x 720mm
Shipping Weight	125kg
Adjustments	Thoracic and lumbar drops: magnetic with adjustable resistance and 15mm of rear-hinged drop Head section and armrest: +40° to -35° Extendable footrest: +200mm Height adjustable footrest: +200mm

ITEM CODE
MC5020 Athlegen Apollo Chiropractic Table

Ph: 03 9878 1399 | sales@accesshealth.com.au | www.accesshealth.com.au