

Plinth 2000

Model 502T Traction Table

Overview:

Designed specifically for lumbar, cervical and cervical in sitting traction procedures. Featuring a friction free rolling top for effective application of traction force and supplied with harnesses, masts and a spreader bar. The table is also available as part of a package along with traction machine and flexion stool.




Standard Features

- UK Manufactured
- Lifetime frame & 5 year parts guarantee
- 16 year average lifespan
- Hydraulic or electric operation
- Comfort foam
- Easy operation
- Low running costs
- Hardwearing anti-microbial vinyl
- Safe working load 250kg
- Retractable castors
- Hoist clearance 130mm
- Harness set (cervical & lumbar traction)
- Friction free rolling top
- Mast & spreader bar
- Traction machine table


Overall Dimensions

Length: 198cm (78") approx
Width: 63cm (25") approx
Height: 47cm - 97cm approx
(18" - 38")approx

Upholstery Dimensions

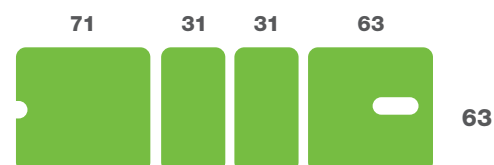

Recommended Options:

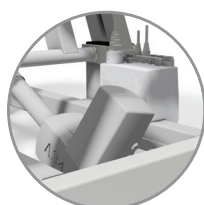

Battery Back-up



Flexion Stool



Traction Machine



Foot Switch



Spare Harness Set

NHS Approved Supplier




t: +44 (0)1449 767887
e: sales@plinth2000.com
www.plinth2000.com



Operational Parameters



Rolling Top Movement: 14cm

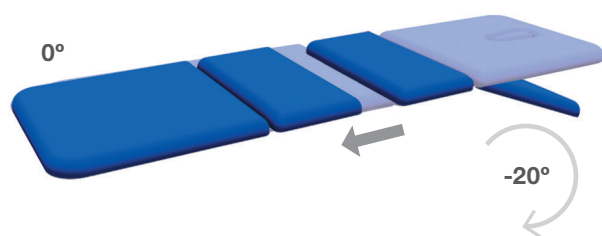

Model 502TE - Electric version

Model 502TH - Hydraulic version

For Your Peace of Mind

- Lifetime frame & 5 year parts guarantee (accidental damage to upholstery excluded)
- Optional extended warranty to 10 years
- Average lifespan of 16 years
- Designed and manufactured in-house in the United Kingdom
- CE marked
- Externally audited to ISO9001:2008 by LRQA

Standards

- IP Rating - IPx4
- EEC Classification - MDD93/42/EEC Annex IX Class 1
- Electrical Standard - IEC601 Part 1 + 60601. P (EN)
- Vinyl Flame Retardant To CRIB 5 (BS7176 Including Foam)
- Foam To BS 5852: Part 2: 1982, Ignition Source 5 [CRIB 5]

Delivery Promise

- Delivery to UK mainland in just 7-10 working days
- Faster service maybe available if required

Upholstery Vinyl

- Flame retardant to CRIB 5
- Patented anti-microbial technology proven to inhibit the growth of E-Coli and MRSA
- Excellent levels of wear and longevity
- Fully washable and impermeable, to aid infection control
- Suitable for cleaning and disinfecting with dilute solutions of Hypochlorite (10,000ppm)
- Please note: Unless specified, beige will be supplied.



Plinth 2000 Ltd.
 Wetheringsett Manor
 Wetheringsett
 Stowmarket
 Suffolk, IP14 5PP

t: +44 (0)1449 767887
 f: +44 (0)1449 766122
 e: sales@plinth2000.com
 www.plinth2000.com

Plinth | 2000

